



SMOBIO

Small Bio, Smart Tool

www.smobio.com

Product Information

ExcelBand™ series

50 bp DNA Ladder, Ready to Use

DM1100 **54 µg / 500 µl**

Size: 500 µl

5 µl per well for 100 gel applications

Storage

Room temperature for 6 months

4°C for 12 months


-20°C for 36 months

Description

The DM1100 ExcelBand™ 50 bp DNA Ladder is a ready-to-use DNA ladder, which is premixed with loading dye for direct gel loading. The DM1100 DNA ladder is composed of 17 individual DNA fragments: 1.5k, 1.2k, 1k, 900, 800, 700, 600, **500**, 450, 400, 350, 300, 250, **200**, 150, 100 and 50 bp derived from a mixture of PCR products and specifically digested plasmid DNA. This product contains two enhanced bands (500 bp and 200 bp) for easy reference. In addition, the low range Orange G tracking dye mimics the migration of a 50 bp dsDNA during electrophoresis, and allows for real time monitoring.

Buffering

10 mM Tris-HCl (pH 8.0), 10 mM Na-EDTA, 0.025% Orange G and 10% Glycerol.

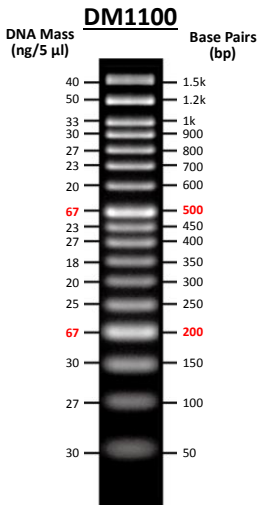

2% Agarose Gel

0.5x TAE Buffer

Note: This DNA Ladder was not designed for precise quantification of DNA mass, but can be used for approximating the mass of DNA by intensity comparison.




SMOBIO
Small Bio, Smart Tool

ExcelBand™ 50 bp DNA Ladder

DM1100

DNA Mass (ng/5 μ l) Base Pairs (bp)

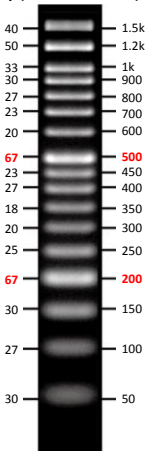

2% Agarose Gel
0.5x TAE Buffer



The ladder is supplied with the gel loading buffer and is ready to use. Do not heat or dilute before loading.

P06-F11-A1

Quality Control

1. 17 defined bands are observed during electrophoresis in agarose gel.
2. No nuclease activity is detected by a direct nuclease activity assay.

Other information

SMOBIO Technology, Inc. claims all warranties with respect to this document, expressed or implied, including but not limited to those of merchantability or fitness for a particular purpose. In no event shall SMOBIO Technology, Inc. be liable, whether in contract, tort, warranty, or under any statute or any other basis for special, incidental, indirect, punitive, multiple or consequential damages in connection with or arising from this document, including but not limited to the use thereof.

Caution: Not intended for human or animal diagnostic or therapeutic uses.

Related Products

DM2100	ExcelBand 100 bp DNA Ladder, 500 μ l
DM2300	ExcelBand 100 bp+3K DNA Ladder, 500 μ l
DM3100	ExcelBand 1 KB (0.25-10 kb) DNA Ladder, 500 μ l
DM5100	ExcelBand Super Range DNA Ladder, (50 bp- 25 kb), 500 μ l
PM2610	ExcelBand Enhanced 3-color High Range Protein Marker, 250 μ l \times 2
PM2800	ExcelBand 3-color Extra Range Protein Marker, 250 μ l \times 2
WM1000	YesBlot Western Marker I, 250 μ l
TP1000	ExcelTaq DNA Polymerase, 500 U
TP1200	ExcelTaq 5 \times PCR Master Dye Mix, 200 RXN
TF1000	SMO-HiFi DNA Polymerase, 100 U
RP1000	ExcelRT Reverse Transcriptase, 20,000 U
DL2000	ExcelDye DNA Loading Dye, Green, 5 ml \times 2
DL4000	ExcelDye DNA Loading Dye, Tri-color, 5 ml \times 2
DL5000	FluoroDye DNA Fluorescent Loading Dye (Green, 6X), 1 ml
DS1000	FluoroStain DNA Fluorescent Staining Dye (Green, 10,000X), 500 μ l
NS1000	FluoroVue Nucleic Acid Gel Stain (10,000X), 500 μ l
VE0100	B-BOX™ Blue Light LED epi-illuminator



B-BOX™ Blue Light LED epi-illuminator

For research use only
2019 ver. 2.1.1

P06-F11-A1